


# Protilátky ke tkáňové transglutamináze

**Anti-tTG IgA, IgG** (atTg) jsou protilátky ke tkáňové transglutamináze. Tkáňová transglutamináza má přímý vztah k patogenezi onemocnění a byla popsána jako vlastní, chemický substrát endomysia. Tkáňová transglutamináza (izoenzym transglutamináza II, TG2 - EC 2.3.2.13), je transferázou, systémový název je protein-glutamin:amin- $\gamma$ -glutamyltransferasa. Je to  $\text{Ca}^{2+}$  dependentní enzym katalyzující deaminaci glutaminu na glutamát, rovněž vede ke vzniku intramolekulární vazby glutaminu na další primární amin, např. lysin, a vede k agregaci glutaminových peptidů. Stanovení protilátek ke tkáňové transglutamináze (atTG) má proto rovněž velmi vysokou diagnostickou efektivitu, podobně jako protilátky proti endomysiu (senzitivita 87–97 % a specifita 88–98 %).


**Stanovení atTG** je prováděno klasickou metodou ELISA, což je pro rutinní diagnostiku technika dostupnější než imunofluorescenční průkaz protilátek proti endomysiu (EmA).

Protilátky atTG lze na rozdíl od EmA stanovovat ve třídě IgA i IgG, což má význam pro nemocné se selektivním deficitem IgA. Metoda byla popsána s použitím antigenu morčete, který je použit ve většině starších souprav. Novější soupravy již používají jako antigen tkáňovou transglutaminázu izolovanou z lidských buněk, z lidských erytrocytů, nebo rekombinantní tTG izolovanou na *E. coli*. Referenční hodnoty se liší u jednotlivých souprav, většinou je pro IgA protilátky uváděna horní hranice normy 10–15 IU/l, některé soupravy definují i tzv. šedou zónu (*gray-zone*) v rozsahu 10–20 IU/l. Stanovení protilátek atTG s lidským rekombinantním antigenem vykazuje nižší falešnou pozitivitu než metody s antigenem morčete. Nejnovější studie porovnávají protilátky třídy IgA a IgG a POCT metodiky stanovení atTG protilátek.

Stanovení protilátek atTG ve třídě IgA je doporučeno jako základní screeningový test pro diagnostiku celiakie.

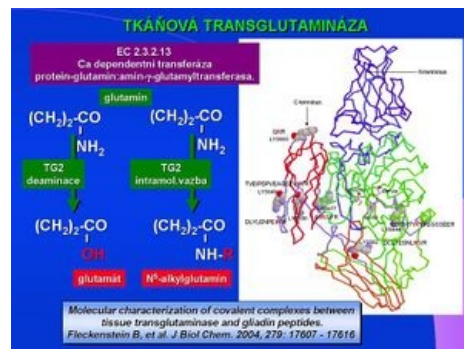

Tkáňová transglutamináza

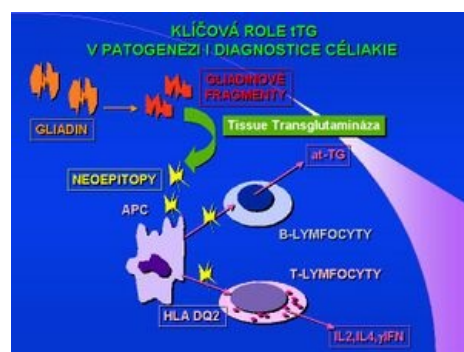

Role tTG v patogenezi a diagnostice celiakie

## Odkazy

### Zdroj

- se svolením autora převzato z KOCNA, Petr. *GastroLab : MiniEncyklopedie laboratorních metod v gastroenterologii* [online]. ©2002. Poslední revize 2011-01-08, [cit. 2011-03-04]. <<http://www1.f1.cuni.cz/~kocna/glab/glency1.htm>>.

### Použitá literatura

- BAVIERA, LC, et al. Celiac disease screening by immunochromatographic visual assays: results of a multicenter study. *J Pediatr Gastroenterol Nutr.* 2007, vol. 45, no. 5, s. 546-50, ISSN 0277-2116 (Print), 1536-4801 (Electronic). PMID: 18030231 (<https://www.ncbi.nlm.nih.gov/pubmed/18030231>).
- SÁRDY, M, et al. Tissue transglutaminase ELISA positivity in autoimmune disease independent of gluten-sensitive disease. *Clin Chim Acta.* 2007, vol. 376, no. 1-2, s. 126-35, ISSN 0009-8981 (Print), 1873-3492 (Electronic). PMID: 16987503 (<https://www.ncbi.nlm.nih.gov/pubmed/16987503>).
- BYRNE, G, et al. Mutagenesis of the catalytic triad of tissue transglutaminase abrogates coeliac disease serum IgA autoantibody binding. *Gut.* 2006, vol. 56, no. 3, s. 336-41, ISSN 0017-5749 (Print), 1468-3288 (Electronic). PMID: 16935926 (<https://www.ncbi.nlm.nih.gov/pubmed/16935926>).
- RAIVIO, T, et al. Self transglutaminase-based rapid coeliac disease antibody detection by a lateral flow method. *Aliment Pharmacol Ther.* 2006, vol. 24, no. 1, s. 147-54, ISSN 0269-2813 (Print), 1365-2036 (Electronic). PMID: 16803613 (<https://www.ncbi.nlm.nih.gov/pubmed/16803613>).
- BARKER, CC, et al. Can tissue transglutaminase antibody titers replace small-bowel biopsy to diagnose celiac disease in select pediatric populations?. *Pediatrics.* 2005, vol. 115, no. 5, s. 1341-6, ISSN 0031-4005 (Print), 1098-4275 (Electronic). PMID: 15867045 (<https://www.ncbi.nlm.nih.gov/pubmed/15867045>).
- MANKÁĚ, A, et al. Tissue transglutaminase antibodies in celiac disease, comparison of an enzyme linked immunosorbent assay and a dot blot assay. *Pathol Biol (Paris).* 2005, vol. 53, no. 4, s. 204-9, ISSN 0369-8114 (Print), 1768-3114 (Electronic). PMID: 15850953 (<https://www.ncbi.nlm.nih.gov/pubmed/15850953>).

- BAUDON, JJ, et al. Diagnosing celiac disease: a comparison of human tissue transglutaminase antibodies with antigliadin and antiendomysium antibodies. *Arch Pediatr Adolesc Med*. 2004, vol. 158, no. 6, s. 584-8, ISSN 1072-4710 (Print), 1538-3628 (Electronic). PMID: 15184223 (<https://www.ncbi.nlm.nih.gov/pubmed/15184223>).
- FERRE-LÓPEZ, S, et al. Immunochromatographic sticks for tissue transglutaminase and antigliadin antibody screening in celiac disease.. *Clin Gastroenterol Hepatol*. 2004, vol. 2, no. 6, s. 480-4, ISSN 1542-3565 (Print), 1542-7714 (Electronic). PMID: 15181616 (<https://www.ncbi.nlm.nih.gov/pubmed/15181616>).
- Tommasini A. - Arch Dis Child. 2004, 15155392. Mass screening for coeliac disease using antihuman transglutaminase antibody assay. *Arch Dis Child*. 2004, vol. 89, no. 6, s. 512-5, ISSN 0003-9888 (Print), 1468-2044 (Electronic). PMID: 15155392 (<https://www.ncbi.nlm.nih.gov/pubmed/15155392>).
- SINCLAIR, D, et al. A comparative study of tissue transglutaminase antibodies and endomysium antibody immunofluorescence in routine clinical laboratory practice. *Ann Clin Biochem*. 2003, vol. 40, s. 411-6, ISSN 0004-5632 (Print), 1758-1001 (Electronic). PMID: 12880544 (<https://www.ncbi.nlm.nih.gov/pubmed/12880544>).
- TONUTTI, E, et al. The role of antitissue transglutaminase assay for the diagnosis and monitoring of coeliac disease: a French-Italian multicentre study. *J Clin Pathol*. 2003, vol. 56, no. 5, s. 389-93, ISSN 0021-9746 (Print), 1472-4146 (Electronic). PMID: 12719462 (<https://www.ncbi.nlm.nih.gov/pubmed/12719462>).
- KOCNA, P, et al. Tissue transglutaminase-serology markers for coeliac disease. *Clin Chem Lab Med*. 2002, vol. 40, no. 5, s. 485-92, ISSN 1434-6621. PMID: 12113293 (<https://www.ncbi.nlm.nih.gov/pubmed/12113293>).
- WOLTERS, V, et al. Human tissue transglutaminase enzyme linked immunosorbent assay outperforms both the guinea pig based tissue transglutaminase assay and anti-endomysium antibodies when screening for coeliac disease. *Eur J Pediatr*. 2002, vol. 161, no. 5, s. 284-7, ISSN 0340-6199 (Print), 1432-1076 (Electronic). PMID: 12012226 (<https://www.ncbi.nlm.nih.gov/pubmed/12012226>).
- HANSSON, T, et al. Recombinant human tissue transglutaminase for diagnosis and follow-up of childhood coeliac disease. *Pediatr Res*. 2002, vol. 51, no. 6, s. 700-5, ISSN 0031-3998 (Print), 1530-0447 (Electronic). PMID: 12032264 (<https://www.ncbi.nlm.nih.gov/pubmed/12032264>).
- CLEMENTE, MG, et al. Antitissue transglutaminase antibodies outside celiac disease. *J Pediatr Gastroenterol Nutr*. 2002, vol. 34, no. 1, s. 31-4, ISSN 0277-2116 (Print), 1536-4801 (Electronic). PMID: 11753161 (<https://www.ncbi.nlm.nih.gov/pubmed/11753161>).
- BARDELLA, MT, et al. Serological markers for coeliac disease: is it time to change?. *Dig Liver Dis*. 2001, vol. 33, no. 5, s. 426-31, ISSN 1590-8658 (Print), 1878-3562 (Electronic). PMID: 11529655 (<https://www.ncbi.nlm.nih.gov/pubmed/11529655>).
- BASSO, D, et al. Role of anti-transglutaminase (anti-tTG), anti-gliadin, and anti-endomysium serum antibodies in diagnosing celiac disease: a comparison of four different commercial kits for anti-tTG determination. *J Clin Lab Anal*. 2001, vol. 15, no. 3, s. 112-5, ISSN 0887-8013 (Print), 1098-2825 (Electronic). PMID: 11344524 (<https://www.ncbi.nlm.nih.gov/pubmed/11344524>).
- CHAN, AW, et al. Tissue transglutaminase enzyme-linked immunosorbent assay as a screening test for celiac disease in pediatric patients. *Pediatrics*. 2001, vol. 107, no. 1, s. E8, ISSN 0031-4005 (Print), 1098-4275 (Electronic). PMID: 11134472 (<https://www.ncbi.nlm.nih.gov/pubmed/11134472>).
- FABIANI, E, et al. The serum IgA class anti-tissue transglutaminase antibodies in the diagnosis and follow up of coeliac disease. Results of an international multi-centre study. International Working Group on Eu-tTG. *Eur J Gastroenterol Hepatol*. 2001, vol. 3, no. 6, s. 659-65, ISSN 0954-691X (Print), 1473-5687 (Electronic). PMID: 11434591 (<https://www.ncbi.nlm.nih.gov/pubmed/11434591>).
- LEON, F, et al. Anti-transglutaminase IgA ELISA: clinical potential and drawbacks in celiac disease diagnosis. *Scand J Gastroenterol*. 2001, vol. 36, no. 8, s. 849-53, ISSN 0036-5521 (Print), 1502-7708 (Electronic). PMID: 11495081 (<https://www.ncbi.nlm.nih.gov/pubmed/11495081>).
- KORDONOURI, O, et al. Autoantibodies to tissue transglutaminase are sensitive serological parameters for detecting silent coeliac disease in patients with Type 1 diabetes mellitus. *Diabet Med*. 2000, vol. 17, no. 6, s. 441-4, ISSN 0742-3071 (Print), 1464-5491 (Electronic). PMID: 10975212 (<https://www.ncbi.nlm.nih.gov/pubmed/10975212>).